


Regular - S

Regular

Wide

- Optimal Interproximal Space
- Combined Diameter of Abutment and Coping
- Smileloc Abutment Diameter

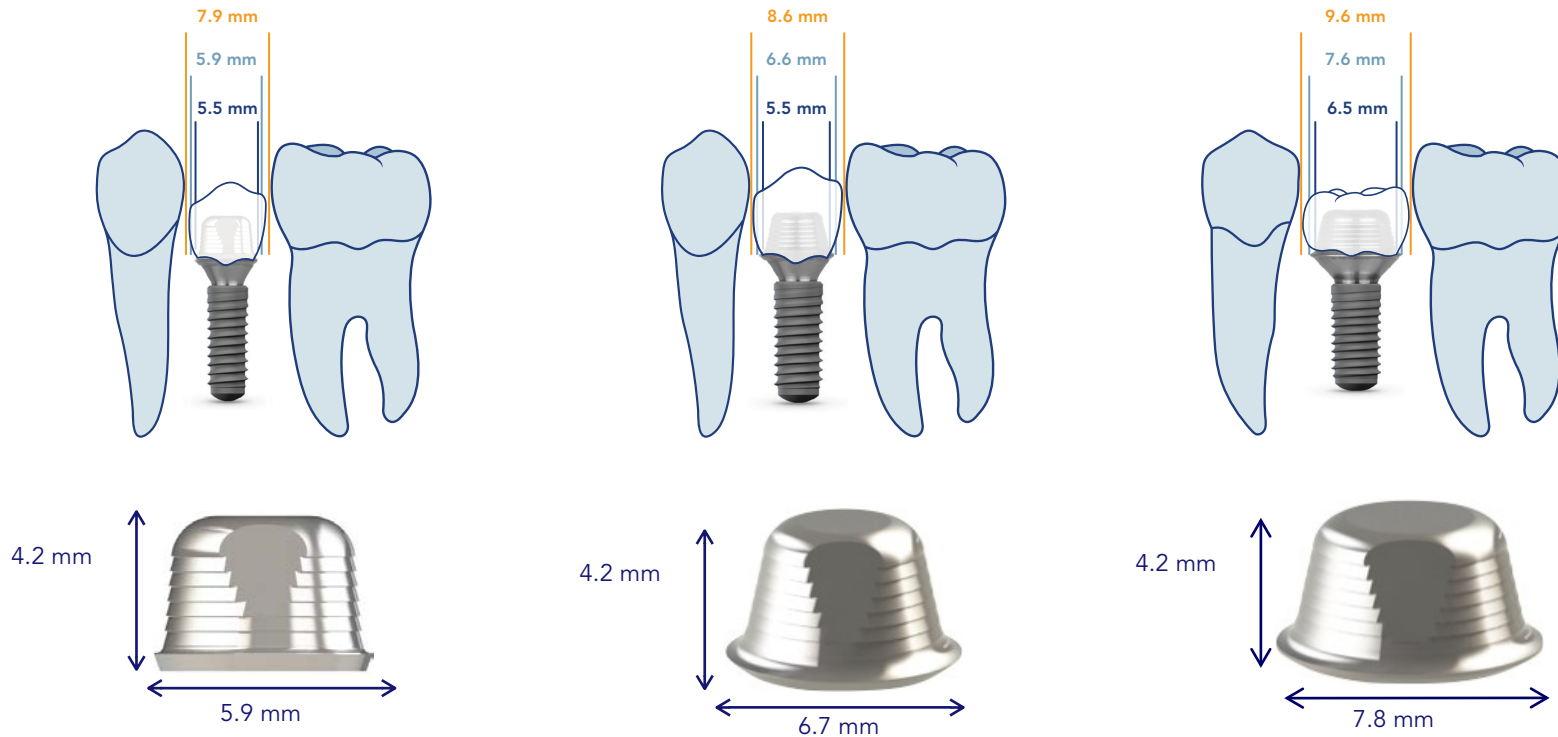

- Minimum interproximal width in the edentulous space:
 - Regular - S Coping: 7.4mm ideally 7.9mm
 - Regular Coping: 8.0mm ideally 8.6mm
 - Wide Coping: 9.0mm ideally 9.6mm
- Implant can be angled up to 10 degrees off the interproximal surfaces of the adjacent teeth.
- A minimum vertical space of 7mm from the top of the implant to the plane of occlusion is required for Zirconia or Emax crowns, unless, a metal occlusal is desired. Reference vertical heights of copings above.